

OASIS WAX FOR STRONG ACIDIC COMPOUNDS

The Oasis WAX (Weak Anion eXchange) SPE material was developed to provide sample preparation for strong acidic compounds. The retention mechanism is mixed mode (both ion-exchange and reversed-phase), which improves retention for strong acidic compounds.

DID YOU KNOW...


You can reduce non specific binding as well as sample loss, when working with therapeutic peptides on μ Elution plates.

Ordering Information

Oasis WAX Sample Extraction Products (Weak Anion Exchange)

Description	Format	Particle Size	Qty.	P/N
Oasis WAX Cartridge	1 cc/10 mg	30 μ m	100/box	186004651
Oasis WAX Cartridge	1 cc/30 mg	30 μ m	100/box	186002489
Oasis WAX Cartridge	3 cc/60 mg	30 μ m	100/box	186002490
Oasis WAX Cartridge	6 cc/150 mg	30 μ m	30/box	186002493
Oasis WAX Cartridge	1 cc/30 mg	60 μ m	100/box	186002491
Oasis WAX Cartridge	3 cc/60 mg	60 μ m	100/box	186002492
Oasis WAX Cartridge	6 cc/500 mg	60 μ m	30/box	186004647
Oasis WAX Plus Cartridge	225 mg	60 μ m	50/box	186003519
Oasis WAX μ Elution Plate	2 mg/96-well	30 μ m	1/pk	186002500
Oasis WAX 96-well Plate	10 mg/96-well	30 μ m	1/pk	186002502
Oasis WAX 96-well Plate	30 mg/96-well	30 μ m	1/pk	186002504
Oasis WAX 96-well Plate	60 mg	30 μ m	1/pk	186003915

Oasis 2 \times 4 Method Development Protocol

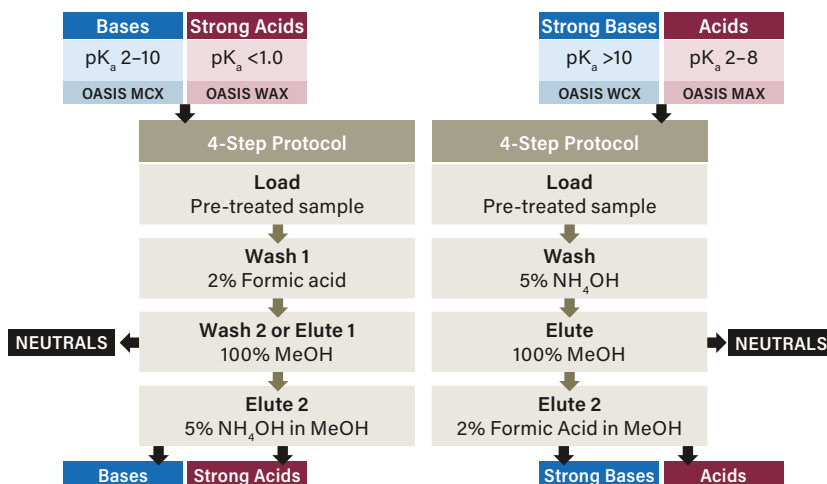

OASIS SORBENT SELECTION TOOLS FOR CONVENIENT METHOD DEVELOPMENT

The Oasis Sorbent Selection Plate and Cartridge Kits enable rapid development of SPE methods for LC-MS analysis. Having all four Oasis Ion-exchange Sorbents (MCX, MAX, WAX, and WCX) in a single plate or a cartridge kit is convenient for scouting the best methods to accomplish efficient isolation of unknown analytes, zwitterionic compounds, or mixtures of analytes with different retention/elution properties.

Ordering Information

Oasis Method Development Kits

Description	Format	Particle Size	P/N
Oasis Sorbent Selection Plate, 3 rows each: MCX, MAX, WCX, and WAX	30 mg/96-well	30 μ m	186003249
Oasis μ Elution Sorbent Selection Plate, 3 rows each: MCX, MAX, WCX, and WAX	2 mg/96-well	30 μ m	186004475
Oasis Sorbent Selection Cartridge Kit, 10 each: MCX, MAX, WCX, and WAX	1 cc/30 mg	30 μ m	186003463
Oasis Sorbent Selection Flangeless Cartridge Kit, 10 each: MCX, MAX, WCX, and WAX	1 cc/10 mg	30 μ m	186006344
Oasis Sorbent Selection Flangeless Cartridge Kit, 10 each: MCX, MAX, WCX, and WAX	1 cc/30 mg	30 μ m	186006345